

Solutions for TGA/STA- EGA coupling

INTRODUCTION

Hyphenated techniques are also known as Evolved Gas Analysis (EGA). They are well established in the field of thermogravimetric analysis (TGA) or simultaneous thermal analysis (STA) with many applications. They have proved particularly interesting for the investigation of gases emitted by the samples enabling their detections, their identifications and in some cases their quantifications.

Hyphenated solutions include MS, FTIR, GCMS and combinations of these.

This technical note describes each of these techniques used in the field of thermal analysis.

1/ Coupling with Mass Spectrometers (MS) – TGA/STA-MS

By coupling a TGA or a STA instrument to a MS, the gases evolved by the sample from the crucible are quickly transferred to the ion source of the mass spectrometer thanks to the pressure difference between the TGA furnace (approx. 1 bar) and the MS (approx. 10⁻⁷ mbar). The gases are not exposed to a temperature lower than 300 °C during the transfer thanks to the use of a heated transfer line.



Figure 1: THEMYS ONE thermal analyzer (on the right) coupled with a mass spectrometer (on the left) from Pfeiffer Vacuum.

After being ionized, the molecules and fragments are then separated according to their mass-to-charge ratio (m/z) thanks to a quadrupole and detected by a SEM or a Faraday detector. The intensity of one or more selected characteristic m/z can be monitored and superimposed to the mass variation signal. Full spectra can also be acquired frequently over the duration of the TGA experiment. Mass spectrometers with amu ranges up to 300 can be employed.

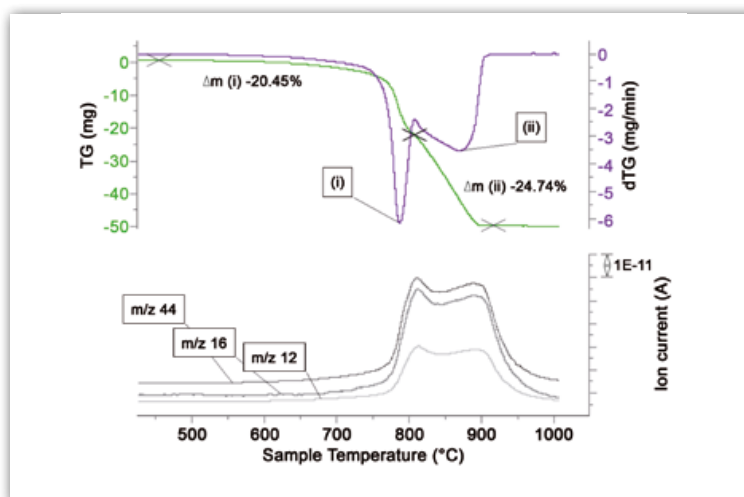


Figure 2 - Two-steps decarbonation of a dolomite (CaMg) (CO₃)₂ sample CO₂ evolution is identified thanks to the ion current intensity of m/z 44, 16 and 12 signals

2/ Coupling with Fourier transformed infra-red spectrometers (FT-IR) – TGA/STA-FTIR

During the coupling of a TGA or STA instrument with a FTIR, the whole carrier and evolved gases from the sample flow through a heated transfer line (up to 300 °C) and up to the FTIR gas cell. Inside the gas cell, chemical groups of the gaseous molecules are identified according to their specific absorption of IR light wavelengths.

The Gram-Schmidt signal - i.e. the integral of the whole spectrum or of a chosen chemical function - can be overlaid with the mass variation. Experimental spectra can be compared with ones from reference materials libraries for identification of evolved substances.

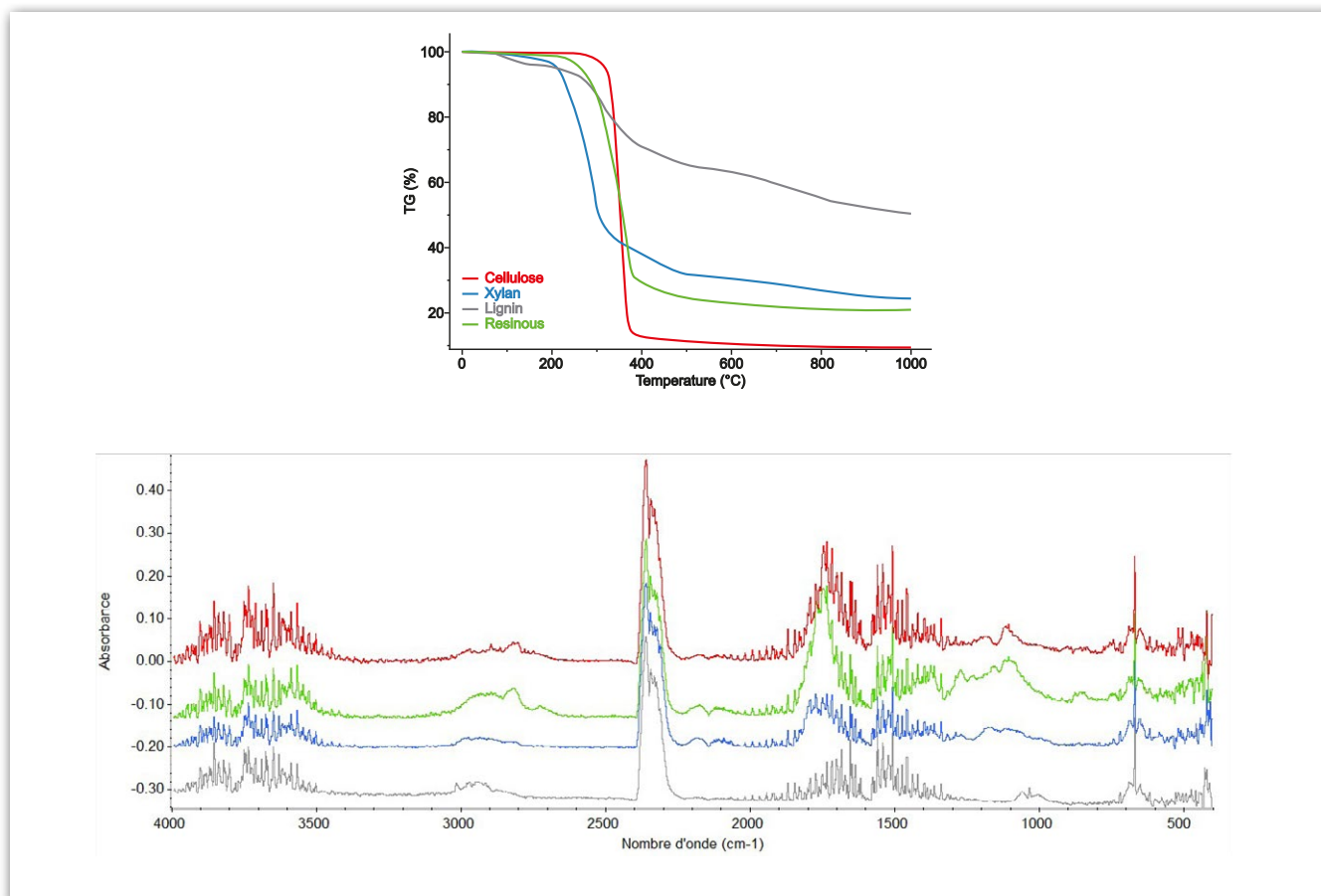


Figure 3 - Biopolymers and wood samples mass losses (up) and corresponding IR spectra at maximum mass loss rates (down). Vibration bands specific to CO₂, H₂O, aldehydes and carboxylic acids could be identified from the cellulose sample IR spectrum.

3/ Coupling with combined gas chromatography and mass spectrometry (GC/MS) - TGA/STA-GC/MS

TGA-GC/MS is becoming increasingly popular, as gas chromatography allows the separation of the evolved species before their detections and identifications by the mass spectrometer. This is particularly interesting when analyzing the thermal decomposition of complex organic substances such as biomass or polymeric materials.

The evolution of the species emitted by the sample from the TGA furnace must be carefully coordinated with the injections on the GC column. This can be accomplished by using an appropriate automated sampling loop that allows the sampling of the whole gas stream from the TGA instrument and their injection on the GC column.

The whole line including the sampling loop has been designed to be heated up to 350 °C to avoid the condensation of the evolved gases from the sample.

Switching between injections and sampling needs to be done in relation to the progress of the sample decomposition process.

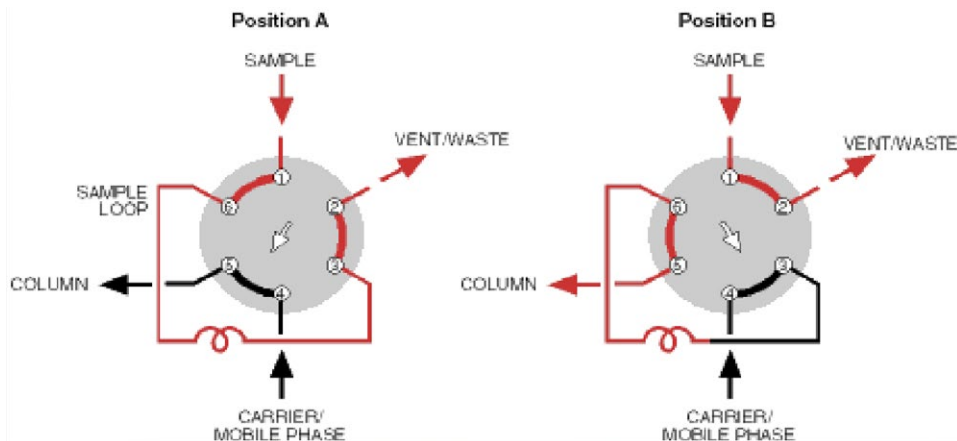


Figure 4: Sample path diagram during sampling phase (position A) and injection into the Autoinjector (position B).

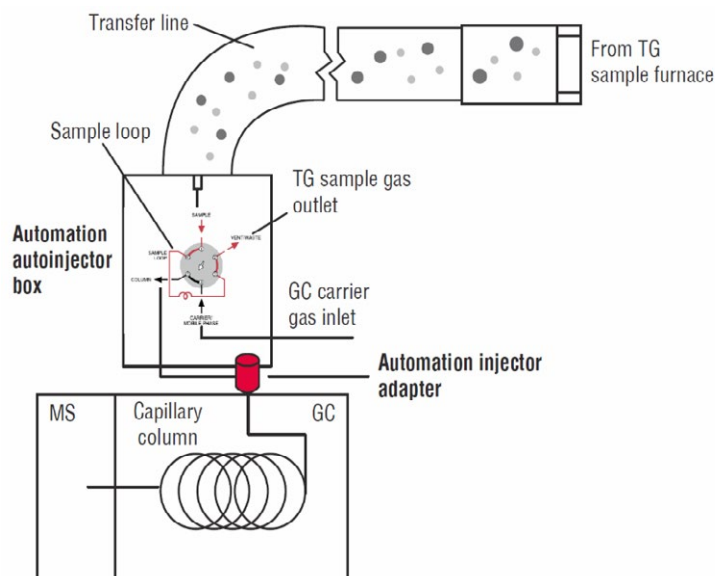


Figure 5: Diagram of the Autoinjector's connection to the GC/MS.

Three modes are available with this Autoinjector:

- Discrete sampling mode: it is meant to achieve a good separation of all the components in the gas mixture but restricts the number of injections that can be analyzed during the TGA run.
- Quasi-continuous sampling mode: the separation cycle is faster and so the gas sampling frequency can be significantly increased. Following the kinetics of production of one specific substance evolved during the TGA run is possible in this mode.
- Timed mode: injections at user-defined times or temperatures corresponding to characteristic TGA events (like mass loss rate maxima).

[Watch video >>](#)

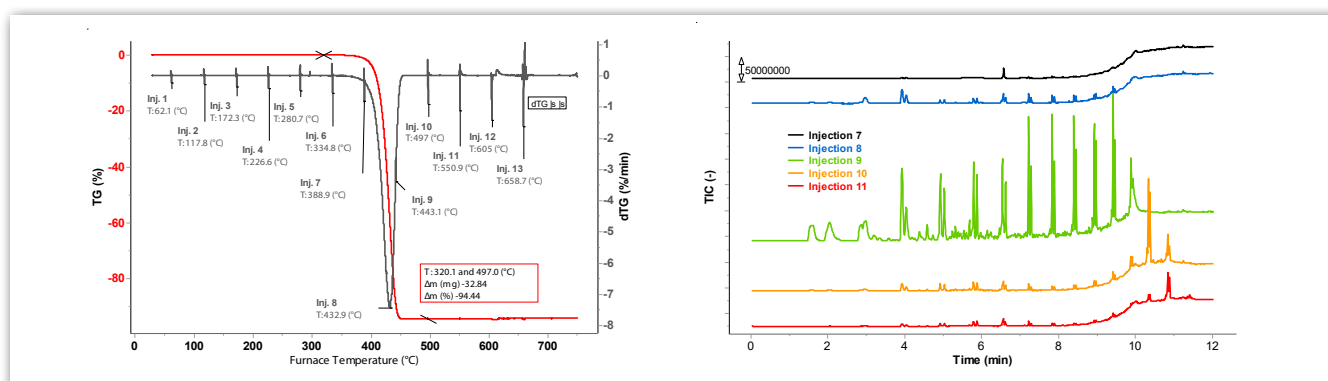


Figure 6 - Mass loss during the pyrolysis of an ethylene and 1-octene copolymer (left) and chromatogram corresponding to injections 7-11 (right). The complex mixture of alkanes and alkenes is separated. Substances like hexadecene (i) and hexadecane (ii) are only detected from injection 10 (497 °C).



Figure 7: THEMYS thermal analyzer (on the right) coupled with a GCMS (on the left) from ThermoFisher.

4. Coupling techniques features & benefits

Technique	Features	Benefits	TGA or STA Models	Temperature range
TGA-MS or STA-MS	Best limit of detection Fast transfer time	Applicable to small mass variations, easy to setup, detection of heavy molecules (with Supersonic System)	CALVET PRO THEMYS ONE / ONE+ THEMYS THEMYS LV	-120 °C up to 2400 °C
TGA-FTIR or STA-FTIR	Signals specific to chemical functions	Applicable to the detection of families of substances in complex gas blends evolved from decomposition of organic samples	CALVET PRO THEMYS ONE / ONE+ THEMYS	-120 °C up to 2400 °C
TGA-GC/MS or STA-GC/MS	Best gaseous molecules separation and identification	Applicable to the identification of substances in complex gas blends evolved from decomposition of organic samples	CALVET PRO* THEMYS ONE / ONE+ THEMYS THEMYS LV*	RT to 2400 °C

*on request


Multiple coupling like TG-GC/MS-FTIR or TG-MS-FTIR are also available for several TGA and STA instruments.

5. Application area of TGA/STA-EGA coupling

TYPICAL APPLICATION AREA	MATERIALS
Gasification, pyrolysis, combustion	Coal, biomass and wastes
Adsorption, desorption, calcination	Catalysts, sorbents, minerals
Dehydration, dehydroxylation, decarbonatation	Cements, minerals, ceramics, pharmaceuticals
Thermal stability, compositional analysis	Pharmaceuticals, polymers, cements, minerals, ceramics
Oxidation	Metals, alloys
Residual solvents	Pharmaceuticals
Polycondensation, curing	Polymers

6. Coupling parts

Coupling part for THEMYS ONE / ONE+

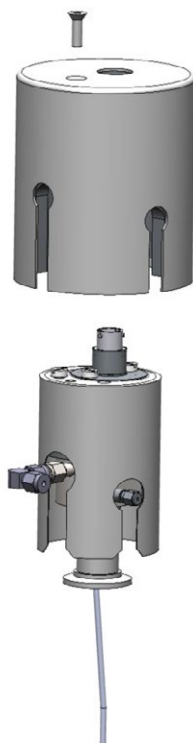

Interface for coupling THEMYS ONE / ONE+ with

FTIR, MS, GC-MS (single coupling)

MS+FTIR, MS+GC-MS (multiple coupling)


Heated up to

300°C

Data acquisition synchronization

Available as an option

Coupling part for THEMYS

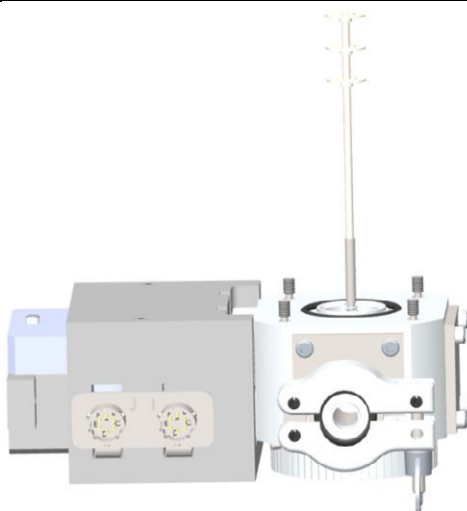

Interface for coupling THEMYS with

FTIR, M S, GC-MS (single coupling)

MS+FTIR, FTIR+GC-MS (multiple coupling)

Heated up to

300°C

Data acquisition synchronization

Available as an option

Coupling part for THEMYS LV

Interface for coupling THEMYS LV with


MS

Heated up to

200°C

Data acquisition synchronization

Available as an option

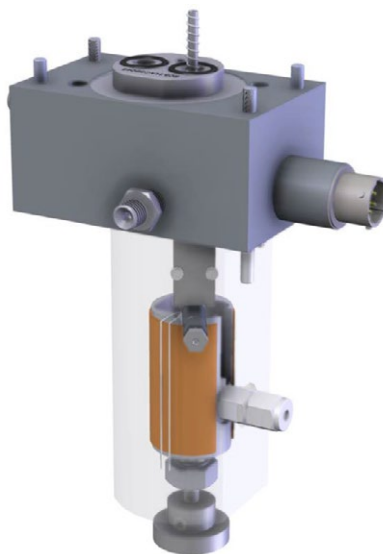Coupling part for CALVET PRO

Interface for coupling CALVET PRO with

FTIR, MS, GC-MS (single coupling)

MS+FTIR, FTIR+GC-MS (multiple coupling)

Heated up to

200°C or 300°C

Data acquisition synchronization

Available as an option